
normavision

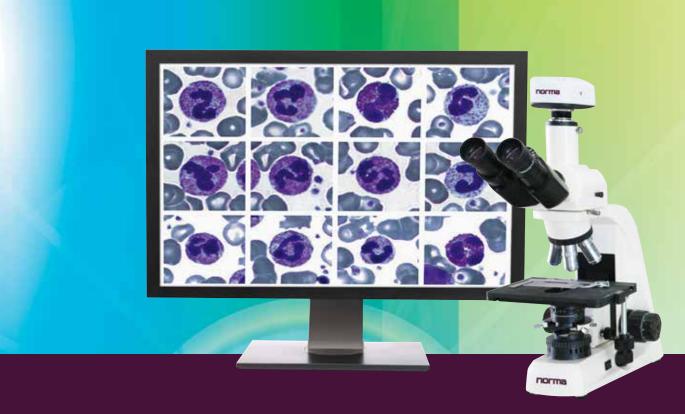


NormaVision M

Analytical blood cell image system


- Digital blood cell morphology for small and medium laboratories
- Help for lab technicians to find and capture blood cell images (WBC, RBC, PLT)
 - Manual identification and automated pre-classification of blood cells
- Validation, database keeping and reports
- Professional development and continuous education

Blood cell analysis process

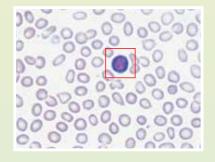

٦

Manual capture of the selected blood cell images (leucocytes, erythrocytes, platelets).

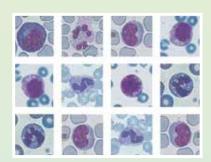

2

After the capturing procedure is completed, the system starts the identification and pre-classification of blood cells.



3

Review and validation of the results by a specialist.



Ч

Integration of the analysis results into the report; consultation.

www.normadiagnostika.com